W

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*The Swedish Back Care System*

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Monteringsanvisning


Assembly instructions

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Instructions de montage

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**Model CN-A, CN-B**

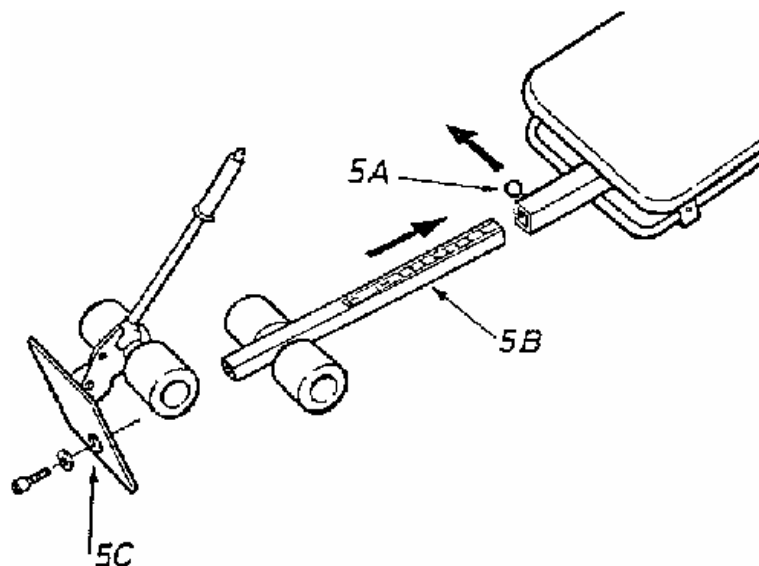Made in China for Mastercare AB, Sweden.

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**English** **ASSEMBLY INSTRUCTIONS - MASTERCARE CN-A, CN-B**

*Important! Always follow the instructions.*

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
**English** Mount the 4 supports with 2 screws in each support (1A). Unfold the pipe framework. Make sure that the steel ropes are fully stretched (1B).

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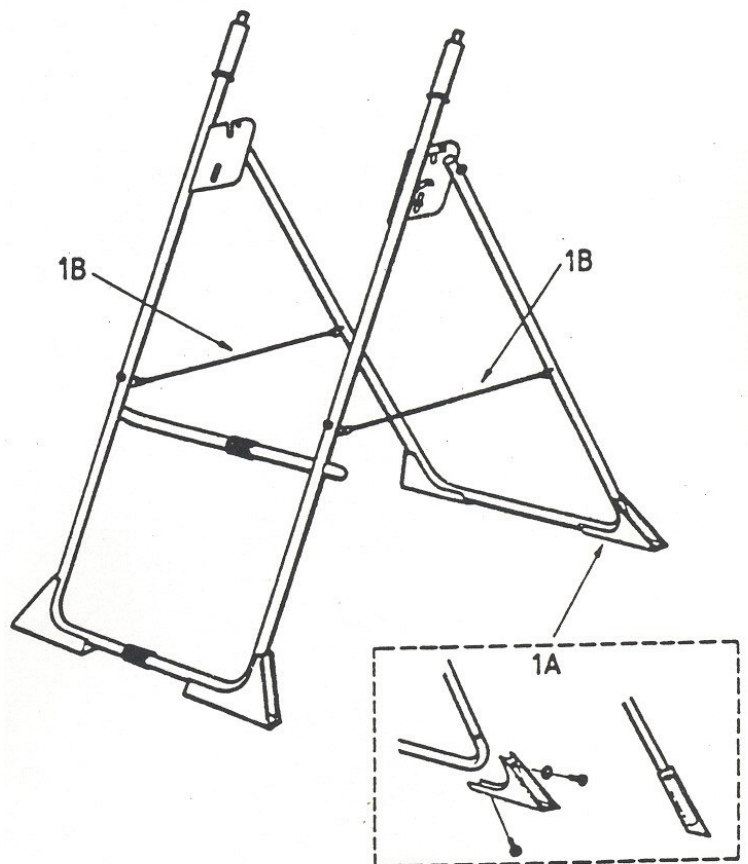English

Attach the suspension devices (2A) one on each side with 3 screws and spring washers on each side.  
Note! The stop pin shall be turned towards the head end of the table (2B).

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English

Put the back rest in its position in the pipe framework (3).

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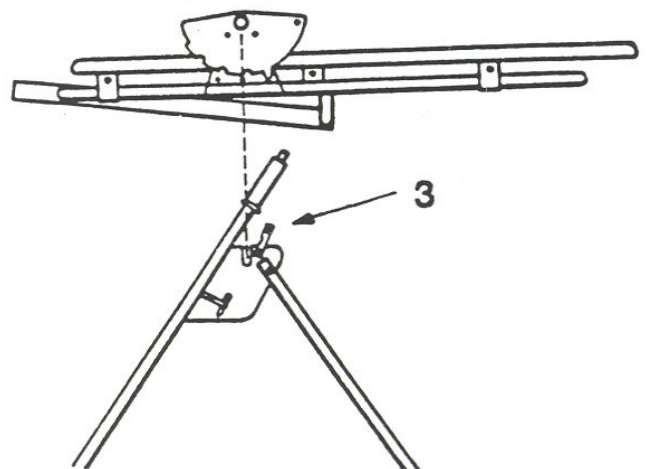
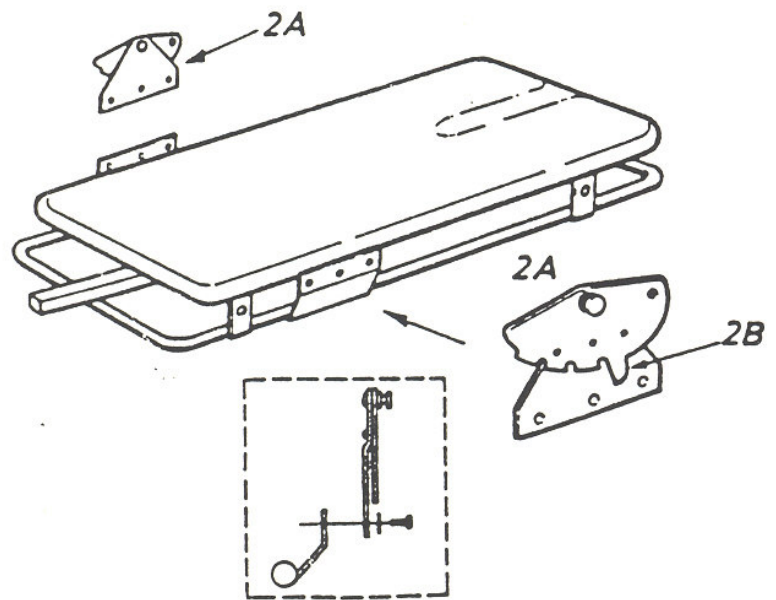
English

Fasten the lock plates (4A). Do not tighten too hard - it might endanger the thread.

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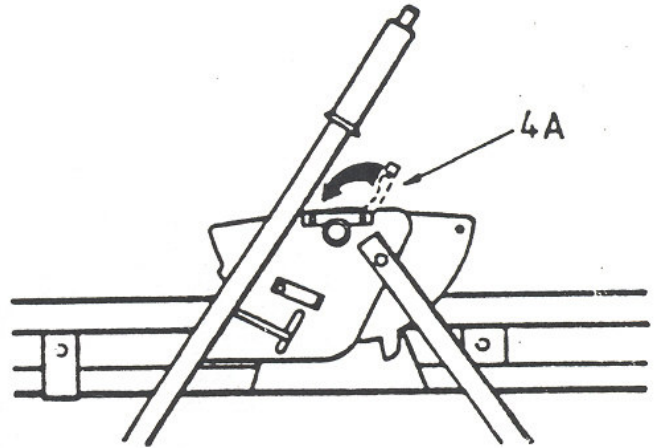
English

Mount the 2 enclosed foot rolls according to this instruction.

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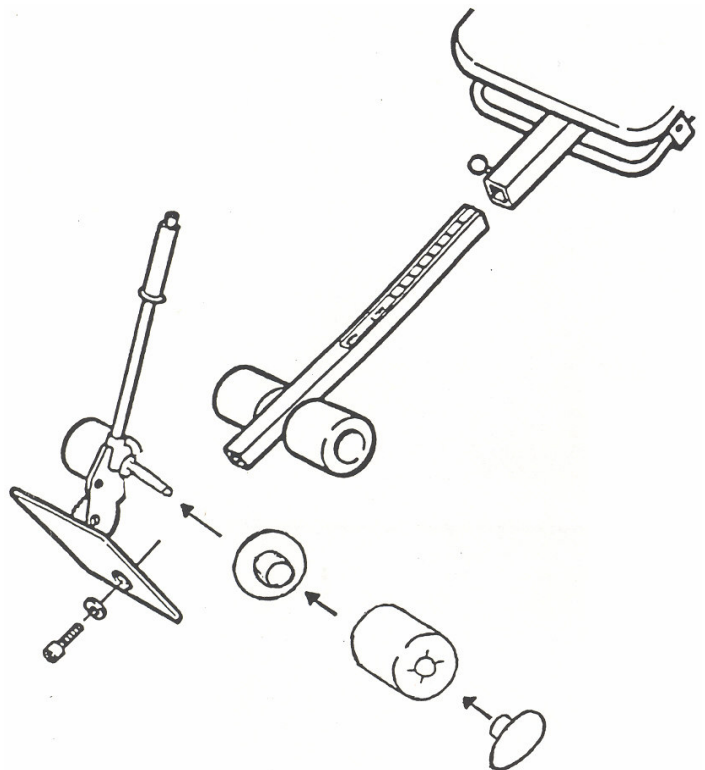
English

Pull the black knob out (5A) and insert the square sectioned pipe (5B).  
Note! The length scale must point upwards. Secure the foot plate by screwing it on (5C).

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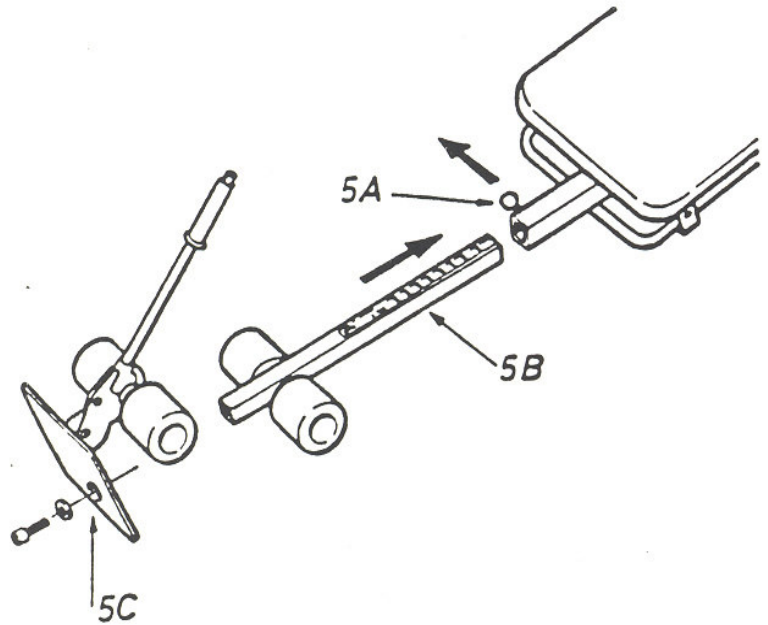
**IMPORTANT!**

Always check and make sure that all screws and nuts are lodging firmly. MASTERCARE must be used on a horizontal base only. Read the instructions carefully before you start using your MASTERCARE.

**PS**

If there is something that you don't understand, please contact your dealer for further information and help.

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